



**Main Features**

- Narrow Band Imaging combined with high-resolution HDTV imaging\* enhances observation of the mucosa and capillaries.  
\* Available only in combination with the CV-180, CLV-180 (CLV-160), and OEV191H.
- Close Focus enables an enlarged, close-up image to be obtained simply by moving the scope tip up to 2 mm from the site.
- Slim design with a diameter of 11.7 mm at the distal end and 11.8 mm at the insertion tube, incorporating a channel with 3.2 mm inner diameter for treatment and suction.
- Highly acclaimed innoflex Variable Stiffness capability enables insertion tube flexibility to be adjusted according to the patient's anatomy.
- Wide 140° degree field of view and extensive four way angulation (180° degree up/down and 160° degree right/left) for observation through all the contours of the colon.
- Auxiliary water jet enables mucus and debris in the colon to be washed away at a touch of a button on the scope or on the footswitch unit when the optional flushing pump is connected.
- Ergonomically designed grip enhances scope maneuverability while easy-to-access controls and user programmable switches improve operability.

- Compatible with the CV-180.
- Scope ID function stores individual scope information in the built-in memory chip and displays it on the monitor, facilitating endoscopy suite management.




Start-up display



**Specifications**

Optical System	Field of view	140°
	Direction of view	Forward viewing
	Depth of field	2 to 100 mm
Distal End	Outer diameter	11.7 mm
	Insertion Tube	Outer diameter
Bending Section	Angulation range	Up 180°, Down 180°, Right 160°, Left 160°
Working Length		L:1680 mm, I:1330 mm
Total Length		L:2005 mm, I:1655 mm
Instrument Channel	Inner diameter	3.2 mm
	Minimum visible distance	3 mm from the distal end
	Endotherapy accessory entrance/exit position in field of view	

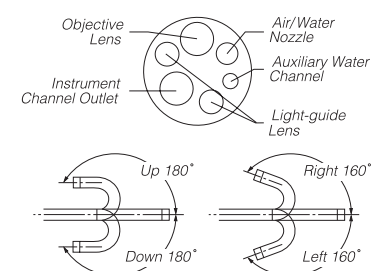

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.



**OLYMPUS MEDICAL SYSTEMS CORP.**  
Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan  
**OLYMPUS MEDICAL SYSTEMS EUROPA GMBH**  
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany  
**OLYMPUS AMERICA INC.**  
3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.  
**OLYMPUS LATIN AMERICA, INC.**  
5301 Blue Lagoon Drive, Suite 290 Miami, Florida 33126-2097, U.S.A.  
**KEYMED LTD.**  
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

**OLYMPUS SINGAPORE PTE LTD.**  
491B River Valley Road #12-01/04 Valley Point Office Tower, Singapore 248373  
**OLYMPUS HONG KONG AND CHINA LIMITED**  
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong  
**OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.**  
Room 1202, NCI Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China  
**OLYMPUS MOSCOW LIMITED LIABILITY COMPANY**  
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia  
**OLYMPUS AUSTRALIA PTY LTD**  
31 Gilby Road, Mount Waverley, VIC., 3149, Australia

www.olympus.com



Your Vision, Our Future

EVIS EXERA II COLONOVIDEOSCOPE  
**OLYMPUS PCF TYPE H180AL/I**



*HDTV, Narrow Band Imaging,  
and outstanding maneuverability  
in an extra-slim colonoscope*

